




# Portable Atomized Ultrasonophoresis Device

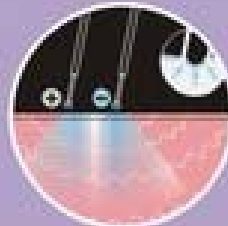

## Instruction Manual

◇ Atomizer ◇ Ultrasonic ◇ Ionic ◇ Photonic

MJ-326

*Before operating the unit, please read this manual completely and retain it for future reference.*

◇ Photonic ◇ Ionic

◇ Atomizer ◇ Ultrasonic

1



**Nano-dynamic Photon**

2



**Iontophoresis Sticks**



**Atomizer**

3



**Ultrasonic Waves**

4

### **01. INTRODUCTION**

When your metabolism is not functioning properly the skin loses moisture and dries out. This creates lines and wrinkles and a loss of elasticity attracting freckles. By using the **Ultrasonophoresis Facial Device** regularly you activate your metabolism which enhances your skin tone and elasticity without using chemicals on your skin or injections. Your blood circulation improves as well as nutrition and vitamin C absorption. Thus creating beautiful youthful skin.

Our human skin consists of three layers; Cuticle, Derma, Subcutaneous Tissue. In the Derma while new cells are produced on a daily basis waste materials also build up. The device removes the build up of waste materials by the process of exfoliation, loosening facial muscles and elasticizing the skin. The overall result is reversing the aging process for your skin. It is a tested and effective way to create skin that has a healthy and radiant glow as well as benefiting from the reduction of lines and wrinkles. MJ-325's settings create visible anti-aging benefits by creating healthier younger looking skin.

## **02. EFFICIENCY**

### **◆ Ultrasonic Waves**

Experience the ultimate in change. Ultrasonic Vibrations (18,000 - 30,000 vibrations per minute.) gives you visible Anti-Aging benefits. The appearance of lines and wrinkles is reduced. Evenness of skin color is improved as age spots and discoloration disappear.

Aged skin will become supple. Since ultrasonic waves stimulate skin muscles with vertical and centrifugal movements, ultrasonic reactivates the aged skin cells, making the skin supple and healthy, making a beautiful face surface. Enhance moisturizer absorption to stimulate cells and revive skin fibers by one million supersonic vibration frequencies per second.

Six kinds of ultrasonic wave forms designed in beauty skin cooperate with mutual changes to improve skin type. Six kinds of ultrasonic wave forms designed in beauty skin, cooperate with mutual changes to improves skin type.

## *Instruction Manual*

---

### ◆ Photonic Waves

To repair cell and skin, red light revive, blue light kill bacteria and yellow light lighten spot.

**Red Photon wave:** 630 nm (nanometer) wavelength - Suitable for any type of skin, enhance blood circulation, make collagen fiber and elasticity fiber of the deep department reconfigure and resume elasticity at the same time.

**Blue Photon wave:** 470nm (nanometer) wavelength - Suitable for sensitive and oil types of skin. Moderate acnitis to achieve getting rid of acne & blemishes. Produce enough collagen , make the cell recombination.

**Yellow Photon wave:** 590nm (nanometer) wavelength - Reduce spot and fleck, restrain melanin.

### ◆ Iontophoresis Sticks


The electrical stimulation with polarity reversal creates temporary openings in the skin's epidermal-dermal barrier, enabling the passage of substances into the dermis for enhancing the penetration of molecules into the skin. With iontophoresis, we enable deeper and more effective penetration of skincare products into the dermis, where it acts on the fibroblasts to stimulate new collagen and elastin productions.

### ◆ Atomizer

An ingenious and convenient design for portable skincare usage. Lotions moisturizers can be poured into the device with built-in an atomizer and directly spray on the face. Then use ultrasonic or iontophoresis to enhance lotions and moisturizer absorption and make your cosmeceutical 100% effect.

### **03. STRUCTURE OF DEVICE**

- 1). Atomizer
- 2). Ultrasonic Metal Probe
- 3). Nano-dynamic Photon
  - red light - red light pulse - blue light - blue light pulse - yellow light - yellow light pulse - red/blue light - red/blue light pulse - red/yellow light - red/yellow light pulse - blue/yellow light - blue/yellow light pulse - red/blue/yellow light - red/blue/yellow light pulse - red/blue/yellow light circle changeable - off.
- 4). Iontophoresis Sticks
- 5). Power Output Indicator LED: Low – Medium – High
- 6). Continuation & Pulse Indicator LED
- 7). Ionic Indicator LED (Red / Green)
- 8). Power On/Off button
- 9). Continuation & Pulse mode button
- 10). Ionic On/Off button
- 11). Photonic On/Off button

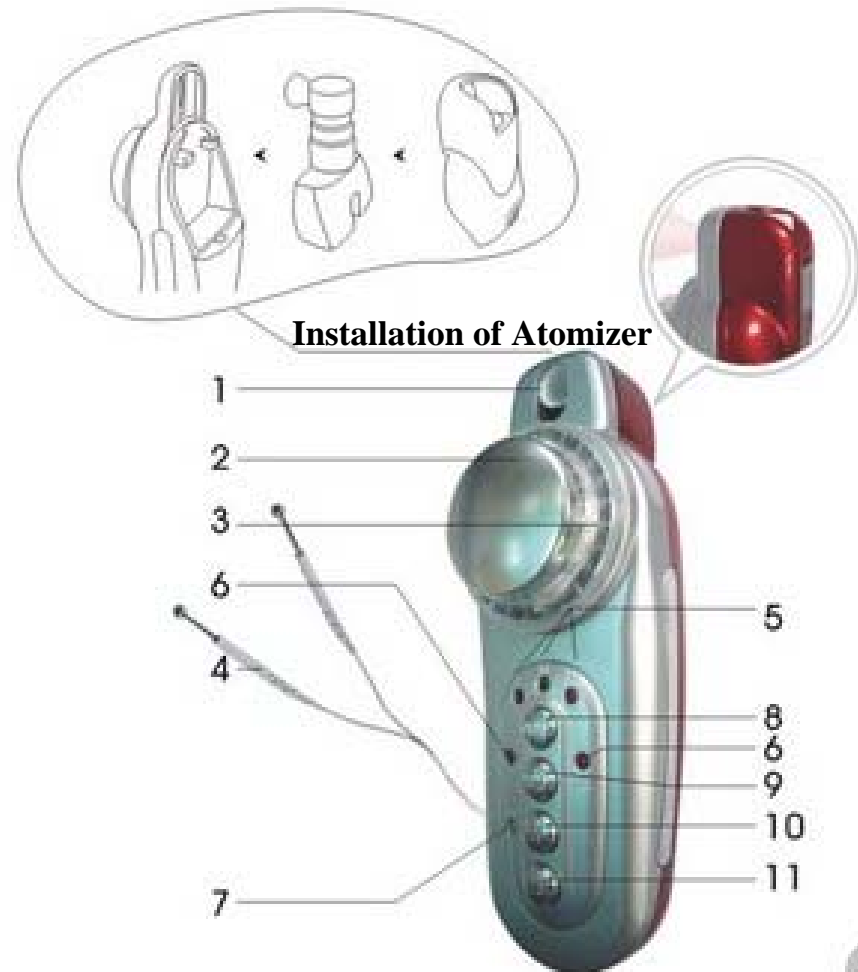

## *Instruction Manual*

### **04. OPERATION OF CONTROL PANEL**

#### ***Power On / Off Button:***

- ◎ This turns on or off the Ultrasonic function, when the Ultrasonic at ON mode, to adjust its weak, medium, strong effect
- ◎ Each press on the button will be accompanied with the correspondent LED light shown.
- ◎ ON → L Speed → M Speed → H Speed → OFF

#### ***Ionic On/Off Button:***

- ◎ When turn on the Ultrasonic Function, the Positive Ionic Function will be on and accompany with the correspondent red LED lighted.  
Notice: Negative Ionic Function is unable to switch on while the Ultrasonic Function is used.
- ◎ Negative Ionic Function can be only turned on when the Ultrasonic Function is off. The correspondent green LED lighted.


#### ***Mode Selection Button:***

Switch between Continuous and Pause mode.

#### ***Photonic On/Off Button:***

- 5 Press the button to switch on Photonic function. Each press to change mode of light.

***Automatic timer setting for 10 minutes, each 10 minutes operating accompanied with bee sound signal.***

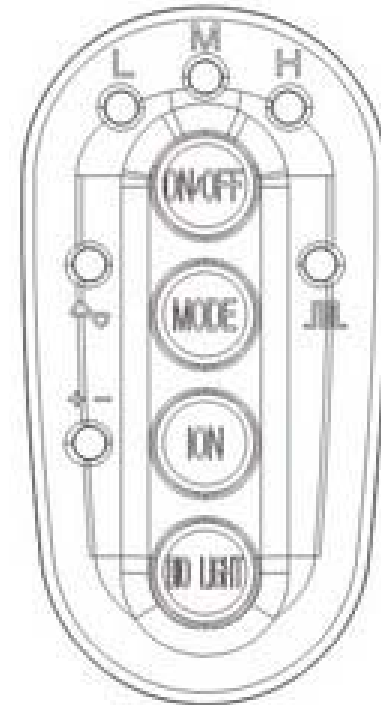

### 05. USE STEPS

(Note: Always apply sufficient skincare products on the face before the treatment. Should feel skincare products are insufficient, please reapply during the treatment.)

#### ◆ Ultrasonic & Photonic Functions:

(Notice: Ultrasonic and Photonic Functions can be used cooperatively or alternatively.)

Step 1 : The cosmetic must be completely wiped off from the face. Clean the face by the Cleansing Gel or Facial Wash.

Step 2 : For further cleanness, use Exfoliating Cream to get rid of dead cells and exfoliate your skin.

Step 3 : Apply softener lotion and softly wipe face, then as suggested to use **Ultrasonic Function** for depth cleaning and removing old dead cells and waste materials; balancing the pH of the skin, constricting skin pores and moisturizing skin.

Step 4 : In order to tenderize and smooth your skin and reviving to the young and balanced skin, please well apply the Firming Gel or Repair Cream on your skin softly.

Step 5 : Then use **Photonic Function** and select continue/pulse and strength of energy functions according to using part to activate your cell and enhance cell metabolism energy.

Step 6 : After the treatment, since cells of skin are actively, please spray the Aliment Cream or Softener Lotion on to moisturize the skin for the skin relaxed and stabilized. Finally, please to wipe the skin care products, and unplug the unit when all treatment completed.

(Notice: By means of ultrasonic vibrations and colored photonic wave accelerate blood and lymph circulation, dispel whole body top layer and deep fat, discharge accumulated toxin and tighten up the skin progressively through strengthening in accordance with difference skincare products.)

## *Instruction Manual*

---

### ◆Iontophoresis Sticks

Step 1. The cosmetic must be completely wiped off from the face. Clean the face by the Cleansing Gel or Facial Wash.

Step 2. To Start the Ionic function, please connect the Ionic Sticks with the Main Body and press the Ionic On/Off Button.

Step 3. Please daubing proper eye renovate essence around eyes. Equability daub from inside to outside.

Step 4. From canthus to trail of eye, moving Iontophoresis Sticks nurse touch point on wrinkle. Stop 3-5 seconds every 1cm alternation. According to the different depth of wrinkle, taking 3-5 times recover nurse to get best result.

Step 5. After the nurse, please disconnect the Ionic Sticks and use soft dressing cotton piece to clean the Iontophoresis Sticks nurse touch point and put it into box.

*(Notice: It is suitable to reduce the crow's-foot, dark circle, detoxifying, firming and wrinkle treatment in accordance with difference skincare products.)*

## **06. NOTIFICATION**

### ***Notification Before Using:***

- ⊙ Before using, please fix the main body and make sure it is securely & tightly fixed.
- ⊙ Please clean metal disk before use.
- ⊙ Please collocate with beauty product, such as beauty gel, cream or lotion.

### ***Notification While Using:***

- ⊙ Please softly and slightly touch skin when use, do not push the unit while using on skin.
- ⊙ Please do not use it in humidity place, like bathroom.
- ⊙ Please avoid eyes and lips area while use the device.
- ⊙ Please operate in circle or horizontal massage direction.
- ⊙ Please select Ultrasonic, Photon or Ion-tophoresis functions according to your facial treatment.
- ⊙ In order to achieve better result, please apply a proper amount of the skincare lotion.
- ⊙ Ultrasonic treatment can not exceed 5 seconds in the same part.
- ⊙ Photon treatment can not exceed 4 minutes in the same part.
- ⊙ Iontophoresis treatment can not exceed 5 seconds in the same part.
- ⊙ The full operation should take no longer than 15 minutes.
- ⊙ Please operate the device with slowly movement and do not stay a same part for a long time.
- ⊙ Please consider the skin condition before use. Inflamed and irritable skin can not be operated.
- ⊙ Frequency of use: suggest using 2-3 times a week, please do not overuse.

## 07. NURSE TIPS



### ***Alimentary Absorption:***

Owing to the high-speed and nanometer ultrasonic vibration, the alimentary cream will be divided into particle instantly. Operating with ***Red Photonic wave*** to enhance the absorption.



### ***Acne Amelioration:***

Deep clean skin, helps to retrieve the original healthy skin by removing the old dead cells and other contaminations. Clean and healthy skin starts with cleansing. (Operating with ***Blue Photonic Wave***)



### ***Rejuvenation:***

By means of deep vibrations of ultrasonic, the cellular tissues strongly repeat stretching and shrinking, accelerate the blood circulation, soften the cellulite and excrete waste. (Operating with ***Red/Yellow Photonic Wave***)

## *Instruction Manual*

---



### ***Reduce spot and fleck:***

The ultrasonic wave can enhance the structure of cell, stimulate metabolism and restrain melanin.

(Operating with ***Yellow Photonic Wave***)



### ***Atomizer:***

Lotions moisturizers can be poured into the device with built-in an atomizer and directly spray on the face. Then use ultrasonic or iontophoresis to enhance lotions and moisturizer absorption and make your cosmeceutical 100% effect.



### ***Iontophoresis Sticks:***

Central nursing or caring in parts. enable deeper and more effective penetration of skincare products into the dermis, where it acts on the fibroblasts to stimulate new collagen and elastin productions.

## **08. FORBIDDEN**

Do not use this appliance under the following condition:

- ⊙ The body temperature over 38 degrees.
- ⊙ With acute diseases.
- ⊙ Heart diseases or with cardiograph or electronic medical apparatus.
- ⊙ High blood pressure, brain clog and heavy diseases.
- ⊙ Pregnant woman.
- ⊙ Menses (suggestion)
- ⊙ Social disease.
- ⊙ Hemophilia.
- ⊙ Heart disease.
- ⊙ Body irritability.
- ⊙ Malignant tumor.
- ⊙ Face nerve disease.
- ⊙ Have hormone for a long time or blood vessels enlarge because of liver functional disorder.

*(Note: Close supervision is necessary when this appliance is used by/or near children, invalids or disabled person.)*

---

## **09. SAFETY**

- ⊙ Please do not use or put it at wet place.
- ⊙ If any visible appearance damage, please do not use it again.
- ⊙ When plug in the power line, please keep hands dry.
- ⊙ Please do not place it where the children can attach.
- ⊙ Please do not disassemble by yourself.
- ⊙ Please avoid eyes area.
- ⊙ Please do not use it on the face for a long time.
- ⊙ Please follow the instruction.

## **10. CLEANING AND MAINTENANCE**

- ⊙ Turn off and unplug the power cord while not in use.
- ⊙ Avoid switch on while not in use, which may cause head over heated or damaged.
- ⊙ Please turn off the power and unplug the unit before cleaning and wiping.
- ⊙ Clean head of devices by alcoholic cotton after using.
- ⊙ To prolong the appliance's lifespan, use a soft cloth with a little detergent to wipe away the dust & dirt from the surface of main body.
- ⊙ To reduce the possibility of electric shock, do not let the Ultrasonic/Photon device and main body be splashed or immersed in water.
- ⊙ Please consult with customer's service immediately if of any question.

## **11. SPECIFICATIONS**

- ◎ **Input:** AC100V-240V 50/60Hz
- ◎ **Output:** DC24V 1.2A
- ◎ **Power Consumption:** 15±5 W
- ◎ **Ultrasonic Frequency:** 3MHz.
- ◎ **3 Speed Ultrasonic Waves:** Low-Medium-High.
- ◎ **Functions:** Ultrasonic wave, Iontophoresis Treatment & Photonic wave.
- ◎ **Timer Setting:** 10 minutes auto off.